Link do produktu: <https://factory-outlet.pl/termometr-bezdotykowy-na-podczerwien-infrared-termometr-kwl-f02-p-501.html>

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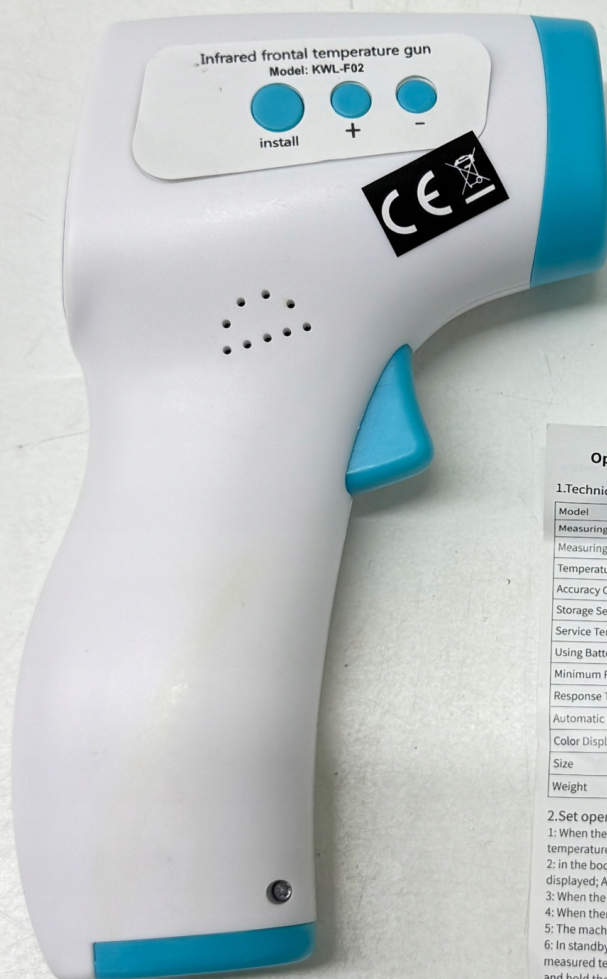
Operation manual of infrared forehead temperature gun

1. Technical parameter

Model	KWL-F02
Measuring Method	Non contact type
Measuring Distance	5cm-8cm
Temperature Measurement Range	Not narrower than 35 °C ~42 °C
Accuracy Of Body Temperature Measurement	35°C~42°C±0.2°C
Storage Service Temperature	-20°C~55°C
Service Temperature	16°C~35°C
Using Batteries	D. C3V (2 No. 7 batteries)
Minimum Reading	0.1°C
Response Time	1S
Automatic Shutdown	<60 seconds
Color Display	Green, red, yellow
Size	150*85*45mm
Weight	105g

2. Set operation

- 1: When the operating range of the ambient temperature exceeds 16 °C ~ 35 °C, the forehead temperature gun will prompt err, or stop providing the temperature value.
- 2: in the body temperature mode, if the measured temperature is lower than 35 °C, "Lo °C" will be displayed; Above 42 °C, the "Hi °C" prompt signal is displayed.
- 3: When the LCD displays "LO", it indicates that the battery is in a low battery state.
- 4: When there is no operation after the test, it will enter the shutdown state within 1 minute.
- 5: The machine has 32 groups of memory for human body and object.
- 6: In standby mode, briefly press the "+" key for 2 seconds to start memorizing and querying the measured temperature. Briefly press the "+" key to increase and the "-" key to decrease. Press and hold the + key and the - key for 5 seconds at the same time to delete all memories
- 7: In standby mode, briefly press "setting key" to switch between "body temperature" and "object temperature" (default body temperature mode when starting up)